



IPL 7800  
COOL PRO

beurer  
beauty

## IPL 7800 COOL PRO IPL hair removal device




Up to 1 mil. light pulses



3.8 cm<sup>2</sup> light area for  
a particularly quick application



5 different intensity levels

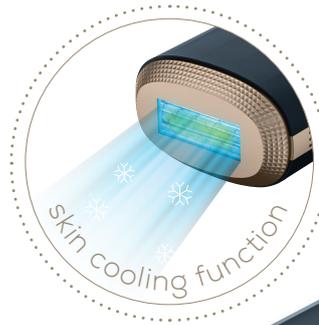

Additional cooling function ensures using the device is as pleasant as possible

State-of-the-art IPL light technology used by dermatologists for professional hair removal. Deactivates the hair roots deep beneath the skin and in doing so decreases the chance of the hair growing back.

3.8 cm<sup>2</sup> light area for a particularly quick application taking around just 23 minutes at the lowest energy level for arms, legs, face, armpits and bikini line

High level of safety with integrated skin contact sensor

Integrated UV filter to protect the skin

With useful „Auto Flash“ mode for convenient application

Up to 50% hair reduction after 3–4 treatments

Cannot be used on dark skin, may be ineffective on white or light-coloured hair.



Mains operation:

Input: 100–240 V~ 50/60 Hz; 1.0 A

Output: 12.0 V === 4.0 A; 48 W

Product measurements (LxWxH): 16.4 x 6.1 x 3.9 mm

Product weight: approx. 280 g

3 year guarantee

Sales unit: 4 / Shipping carton: -

EAN no.: 4211125 57503 8


Item no.: 575.03




## MODE OF ACTION



1 Melanin in the hair absorbs the light.

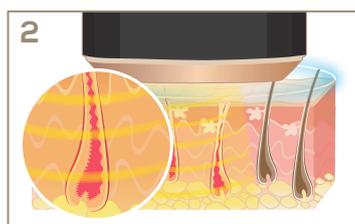

2 The warmth generated deactivates the hair follicle. Additional cooling function ensures using the device is as pleasant as possible.

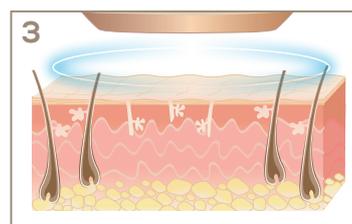

3 The hair falls out and does not continue to grow.

## HAIR AND SKIN COLOUR CHART

		Hair colour					
		Lightest	Light	Medium-Light	Medium	Dark	Darkest
Hair colour	Very Light	X	X	X	X	X	X
	Light	✓	✓	✓	X	X	X
	Medium-Light	✓	✓	✓	✓	X	X
	Medium	✓	✓	✓	✓	✓	X
	Dark	✓	✓	✓	✓	✓	X

## APPLICATION

- ① **Weeks 1 - 2:** use 3 times a week →
- ② **Weeks 3 - 5:** use twice a week →
- ③ **Week 6:** use once or twice a month →

Repeat treatment if required