

EM 70
Wireless

beurer
medical



EM 70 Digital TENS/EMS with remote control unit



TENS
technology

Transcutaneous electrical
nerve stimulation (TENS)



EMS
technology

Electrical muscle stimulation (EMS)



remote
control

Includes remote control unit

PAIN
free

german
engineering



Functions:

Transcutaneous electrical nerve stimulation (TENS),
medicine- and side effect-free pain-relief,
electronic muscle stimulation (EMS) for muscle
training- and regeneration, Massage for relaxation and well-being

4 self-adhesive electrodes including: 1 large electrode (190 x 95 mm),
1 medium electrode (110 x 70 mm) and 2 small electrodes (50 x 50 mm)

Incl. remote control unit

Incl. electrode expansion

19 pre-programmed applications (TENS/EMS/Massage)

Adjustable intensity via device or remote control

Countdown timer

Safety switch-off

Output voltage: max. 90 V_{p-p} (1000 Ohm load)

Output current: max. 90 mA p-p (1000 Ohm load)

Output frequency: 2-290 Hz

Pulse width: 30-300 µs per pulse

Pulse form: dual phase square pulses

"beurer PainAway" app available for free download, App helps you with application:
offers suggestions for electrode placement and explanations of how TENS/EMS work

Device: Rechargeable battery 250 mAh, 3.7 V

Remote Control: 2 x 1.5 AAA batteries (incl.)

Medical product

Product measurements (LxWxH):

Device: 72 x 50 x 14 mm, remote controller: 127 x 49 x 24 mm


Product weight: Device: approx. 38 g, remote controller: approx. 78 g (incl. batteries)

5 year guarantee

Sales unit: 4 / Shipping carton: 4 x 4

EAN no.: 4211125 64821 3

Item no.: 648.21

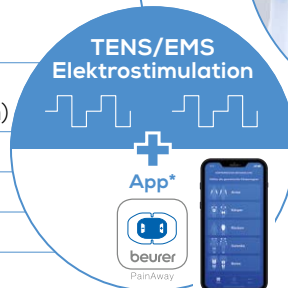

Drug Free
Pain Relief

648.21-EM70_2020-06-15-01_DS_EN_BEU Subject to errors and changes



Remote control included



1. Selection of the body region
2. Selection of the type of pain or muscle building program
3. Setting the intensity



The follow app features await you in the app:

- Helpful explanations about the various Beurer EM 70 Wireless device application programs
- Useful descriptions of the correct electrode positioning on each part of the body
- Instructions on correct usage of the Beurer EM 70
- Additional information concerning the topic of TENS & EMS

The device is controlled exclusively via the separate remote control, so that no connection from the app to the device is required.



The app is available in the following languages:

- German
- English
- French
- Spanish
- Italian
- Turkish
- Russian
- Polish
- Dutch
- Portuguese
- Danish
- Swedish
- Norwegian
- Finnish
- Czech
- Hungarian
- Arabic
- Slovakian
- Chinese



648.21.EM70_2020-06-15_01_DS_EN_BEU Subject to errors and changes

Apple and the Apple logo are trademarks of Apple Inc., registered in the U.S. and other countries. App Store is a service mark of Apple Inc., registered in the U.S. and other countries.

Google Play and the Google Play logo are trademarks of Google LLC. Android is a trademark of Google LLC.



Syncs with
• iOS ≥ 12.0
• Android™ ≥ 8.0



List of supported devices



Drug Free
Pain Relief